Vibrating MassageRoller

Art.nr 373970

NO - DK - SE - GB - DE - FI

Bruk

Enkel å bruke, kun en knapp.

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GB

Simple to use, single button.

Press once to start at low vibration level. M1 will light up.



Press twice from start for medium vibration level. M2 will light up.

Press 3 times from start for high vibration level. M3 will light up.

Press 4 times from start for auto vibration level. M1, M2 and M3 will light up.

Press 5 times from start to switch off. All lights will go out. The device will automatically switch off after 15 minutes of use.

Can be used standing, lying down or sitting.

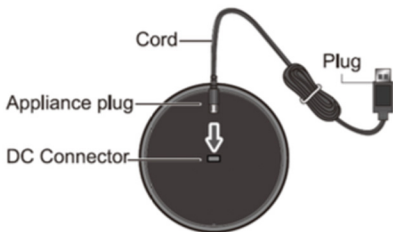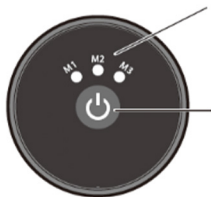Charging & battery

The chargeable lithium battery can be charged using the USB cable supplied. Charge the battery to full capacity for first time use.

When capacity is under 30%, M1 will flash.

When capacity is 30-70%, M1 and M2 will flash.


When capacity is 70-99%, M1, M2 and M3 will flash.

When the battery is fully charged, M1, M2 and M3 will flash.

If the battery does not light after 6 hours of charging, disconnect the product from the charger to avoid danger.

Charging time: approx. 2-2.5 h.

Running time: approx. 2 h.

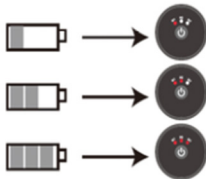Warning

Never leave the product to charge unattended.

Always disconnect the product from the charger when fully charged, before use and before cleaning.

This product is only intended for private use and for use in private households.

It is not watertight, and must not be put in water or any other fluids, left in the rain or in humid conditions.

Keep the product away from children, pets and sharp objects.

Do not use the product or charger if damaged, not working properly or has been dropped on the ground.

Do not use under blankets, duvets, pillows or other textiles.

Any use of training equipment implies a certain risk. Consult your doctor if you are uncertain whether you should use this product.

Troubleshooting

1) Device does not work

- The device is either not charged up sufficiently, or the charger plug may not be inserted properly.

Solution: Charge the device with the USB cable supplied.

- Device not switched on.

Solution: Press the On button.

- The device will automatically switch off after 15 minutes of use.

Solution: Press the On button again.

2) Device gets hot, or cannot be used

- The product has been used repeatedly. Solution: If the device is being charged, disconnect the charger and allow it to cool.

Always disconnect the device from the charger when fully charged.

Cleaning & storage

Clean using a dry or damp cloth.

Store in a dry place, preferably in its box.

IMPORTANT REGARDING SERVICE

In the event of problems of any kind, please contact Mylna Service.

We would like you to contact us before contacting the store so we can offer you the best possible help. Visit our website

www.mylnasport.no - here you will find information about the products, user manuals, a contact form to get in touch with us and possibility to order spare parts. By filling out the contact form you give us the information we need to help you as effectively as possible.

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Mylna Sport AB